


FusionFrame™

Ring Lock System

Quick Connect Struts


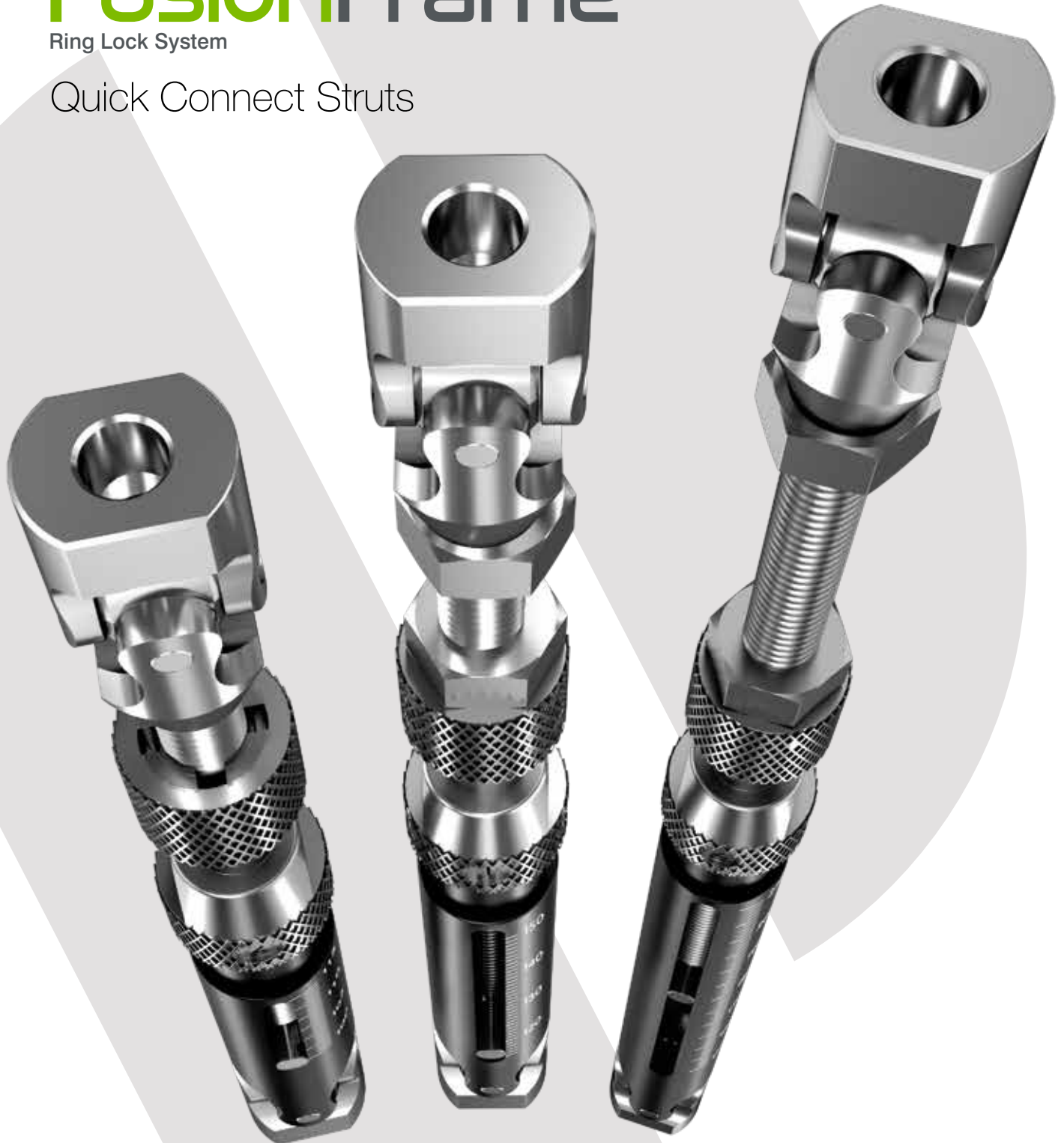


FusionFrame™ Quick Connect Struts

The Quick-Connect Struts provide a mechanically stable and reproducible means for achieving controlled adjustments for the purposes of fracture management and deformity corrections.

- Acute Reductions
- Acute Deformity Corrections
- Gradual Compression and Distraction
- Stabilization of Osteotomies
- Fine-tuned alignment adjustments


Key Features at a Glance

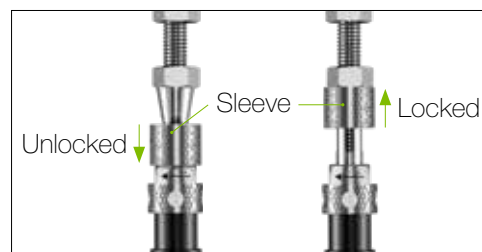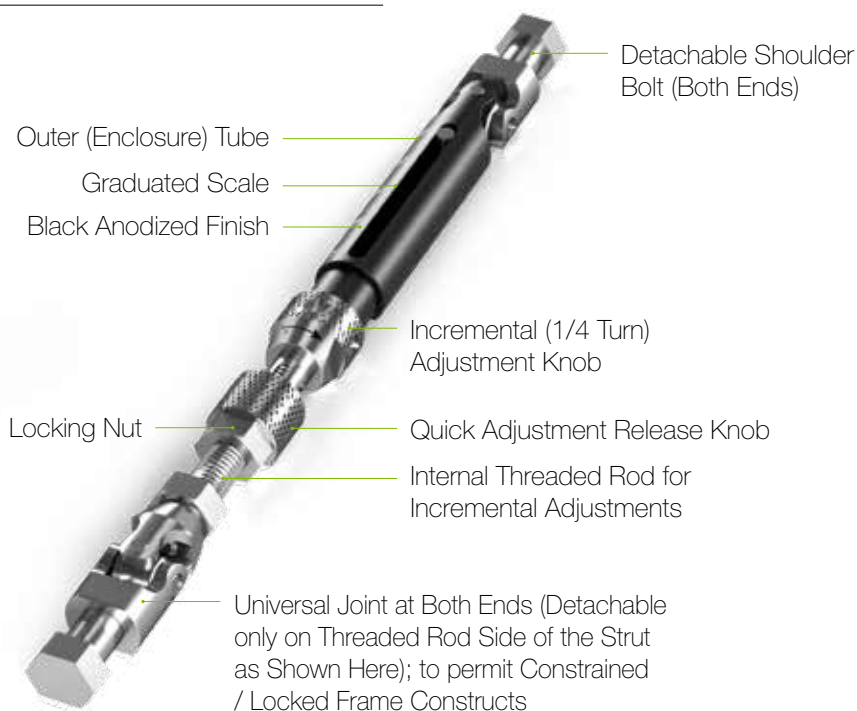

Switching between controlled and quick adjustment is done by pushing / pulling the sleeve over the threaded segment.



Length measuring scale.

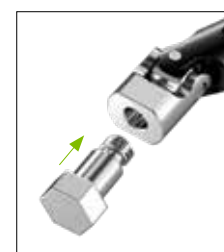

Interface is compatible with Ring holes.

Rapid Length Adjustments

The Struts permit gradual and rapid length adjustments and feature universal joints on both ends to perform acute deformity corrections and assist with fracture reductions.

The universal joint on the threaded end of the Quick Connect Strut is detachable to allow constrained (locked) frame assemblies. A built-in counting knob facilitates the patient's ability to make accurate, one-quarter turns, achieving precisely measured 1/4mm adjustments in the gradual compression or distraction modes.

Strut Size	Fully Collapsed Length (mm)	Fully Distracted Length (mm)	Total Distraction Length (mm)
100mm	100	115	15
116mm	116	152	36
150mm	150	210	60



Distributed by:

innov8ortho, LLC
300 Sylvan Ave, 2nd Floor
Englewood Cliffs, NJ 07632
custsvc@innov8ortho.com

www.glwmed.com

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a surgeon. Rx only.

GLW, Inc., GLW Medical Innovation, GLW Charcot and FusionFrame are trademarks of GLW, Inc.. Innov8ortho is the exclusive distributor of the FusionFrame Circular Fixation System.

© 2023 GLW, Inc. All rights reserved.
Patents: www.glwmed.com/patents/

Literature No: NQQS-SB-V1 DOR: 1020