



creed™
Ortholucet Implants

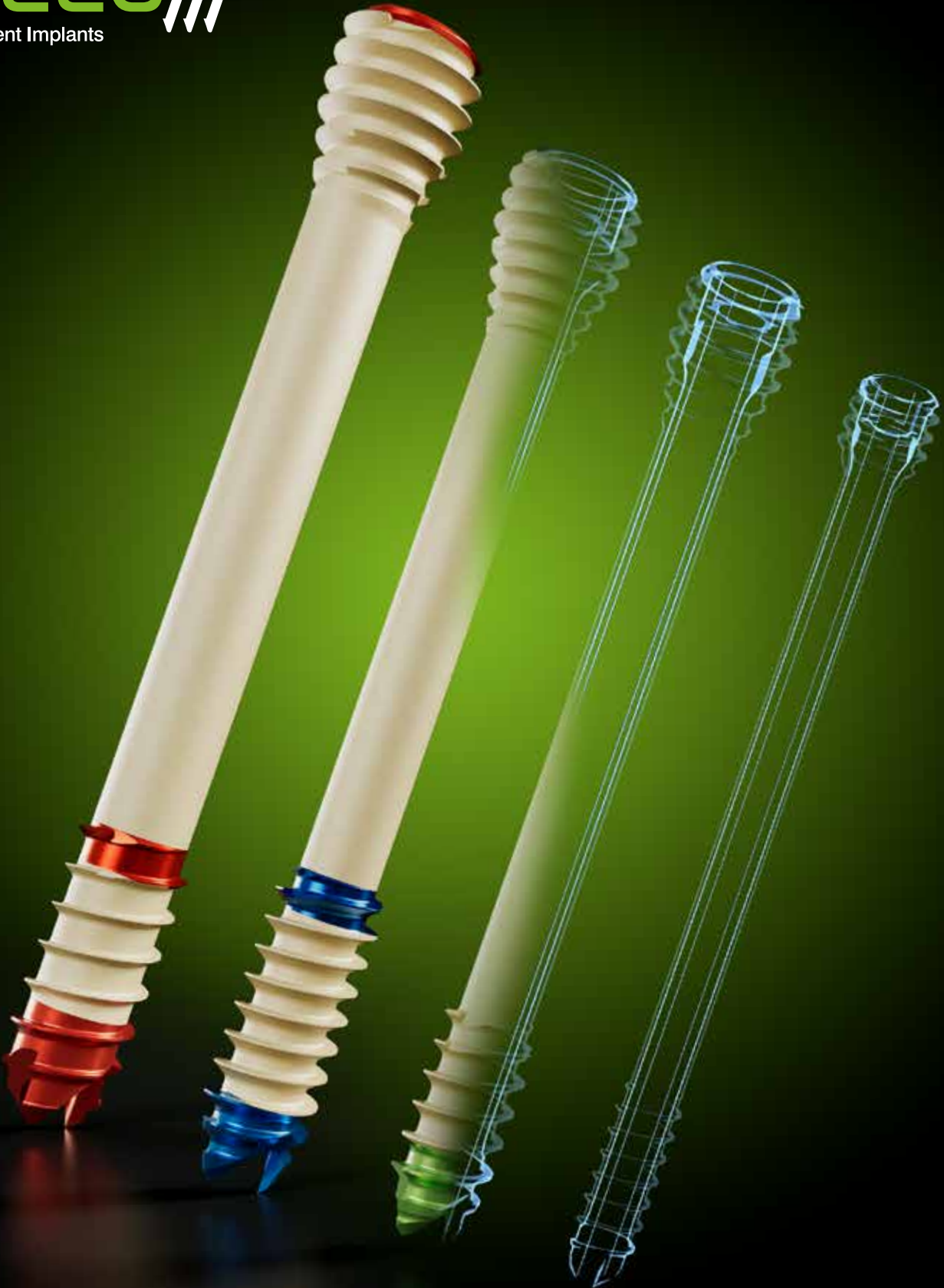


Seeing is **Believing**

GLW Foot & Ankle

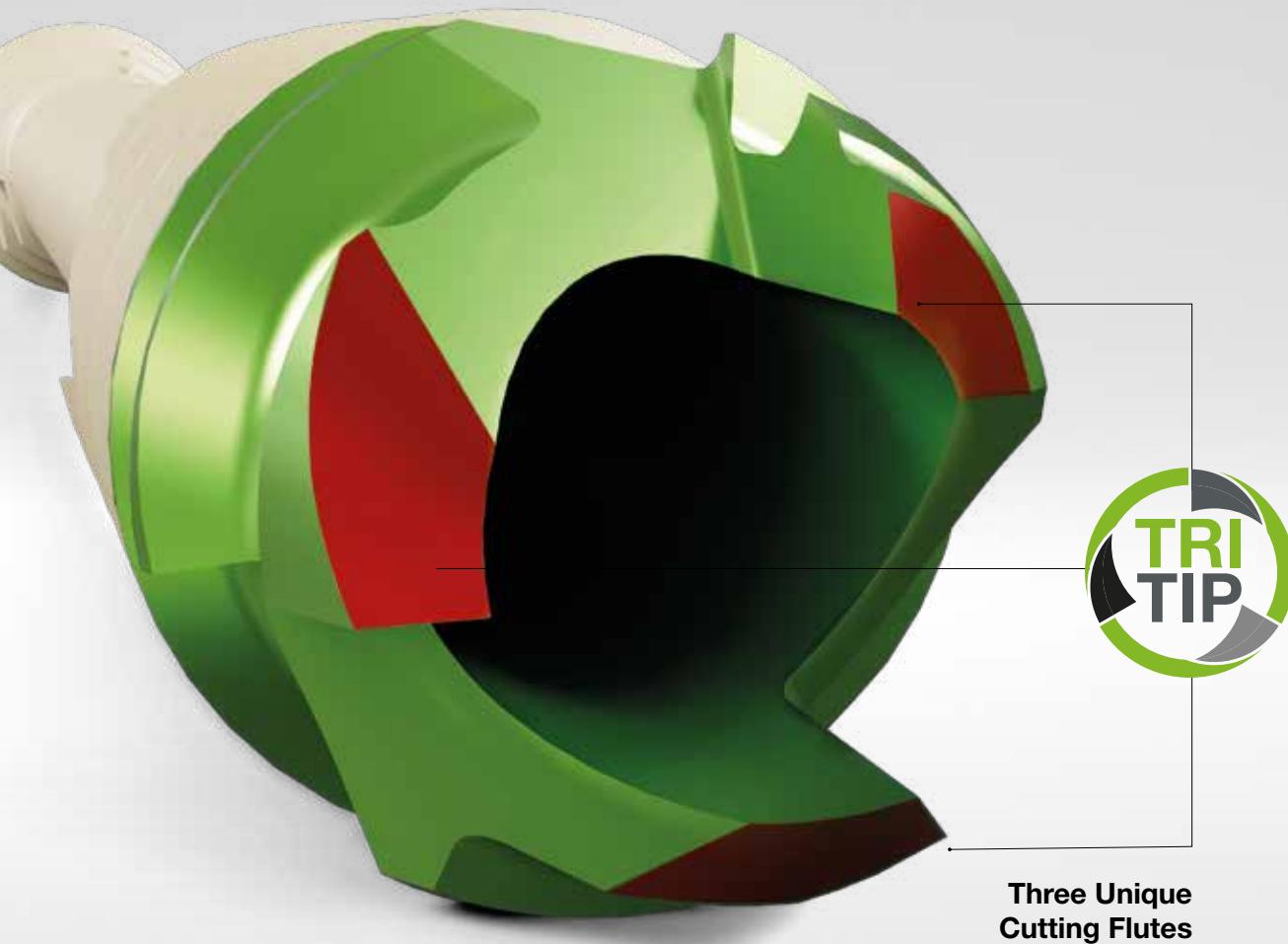

Precision Cutting

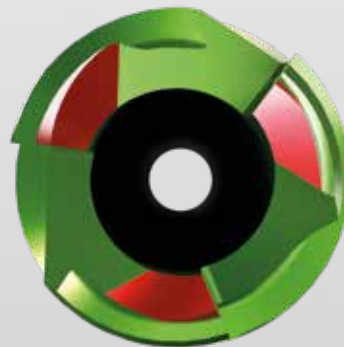

The three aggressive cutting flutes of the screw tip are designed to deliver maximum compression¹ with minimum torque. Superior cutting geometry delivers optimal bone chip displacement, reduces insertion force, and requires no predrilling.²



¹Data on file.

²For Ø7.4mm screw minimized predrilling recommended.


Reverse-Cutting Design
Expedites screw extraction

Forward-Cutting Design
Accelerates screw head countersinking


Ti Core with
Anti-rotation hex provides a strong and stable construct

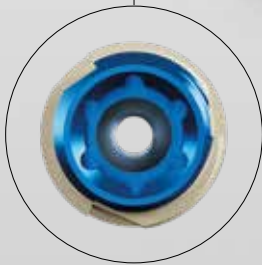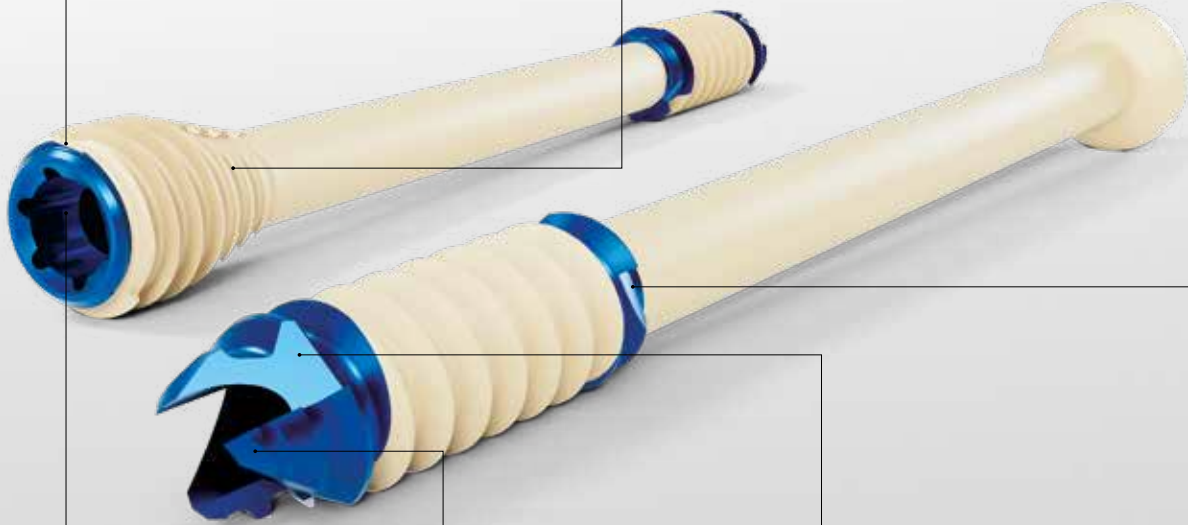

Torx Interface
Promotes a sturdy driver/screw interface

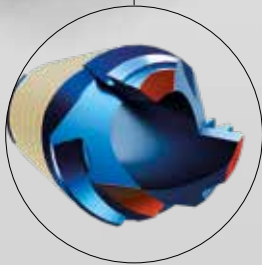


Self-Drilling Tip Design

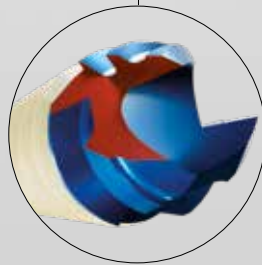

Self-Tapping Tip Design

Streamlines screw insertion

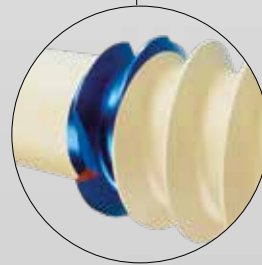

Titanium Reverse-Cutting Feature
Facilitates removal of screw



Stinger Tip


Largest Screw Cannulation

With the largest cannulation and k-wire diameter per screw size on the market, Creed screws are easily inserted using low torque. The larger k-wire, combined with “Stinger Tip” technology, results in a more precise screw insertion trajectory.


Headed	Headless	Ø Screw Size	Ø Screw Cannulation	Ø K-wire Size
		Ø2.5mm	Ø1.4mm	Ø1.3mm
		Ø4.3mm	Ø2.15mm	Ø2.0mm
		Ø5.6mm	Ø2.9mm	Ø2.8mm
		Ø7.4mm	Ø3.3mm	Ø3.2mm

Ortholucen[®] Technology

The patented³ combination of thin-walled titanium, reinforced with a PEEK overmold, provides significant clinical benefits over traditional metal implants by drastically improving the visualization of bones and joint spaces. This novel hybrid material combination provides clear bone visibility during plain radiography as well as less scatter with advanced imaging techniques.

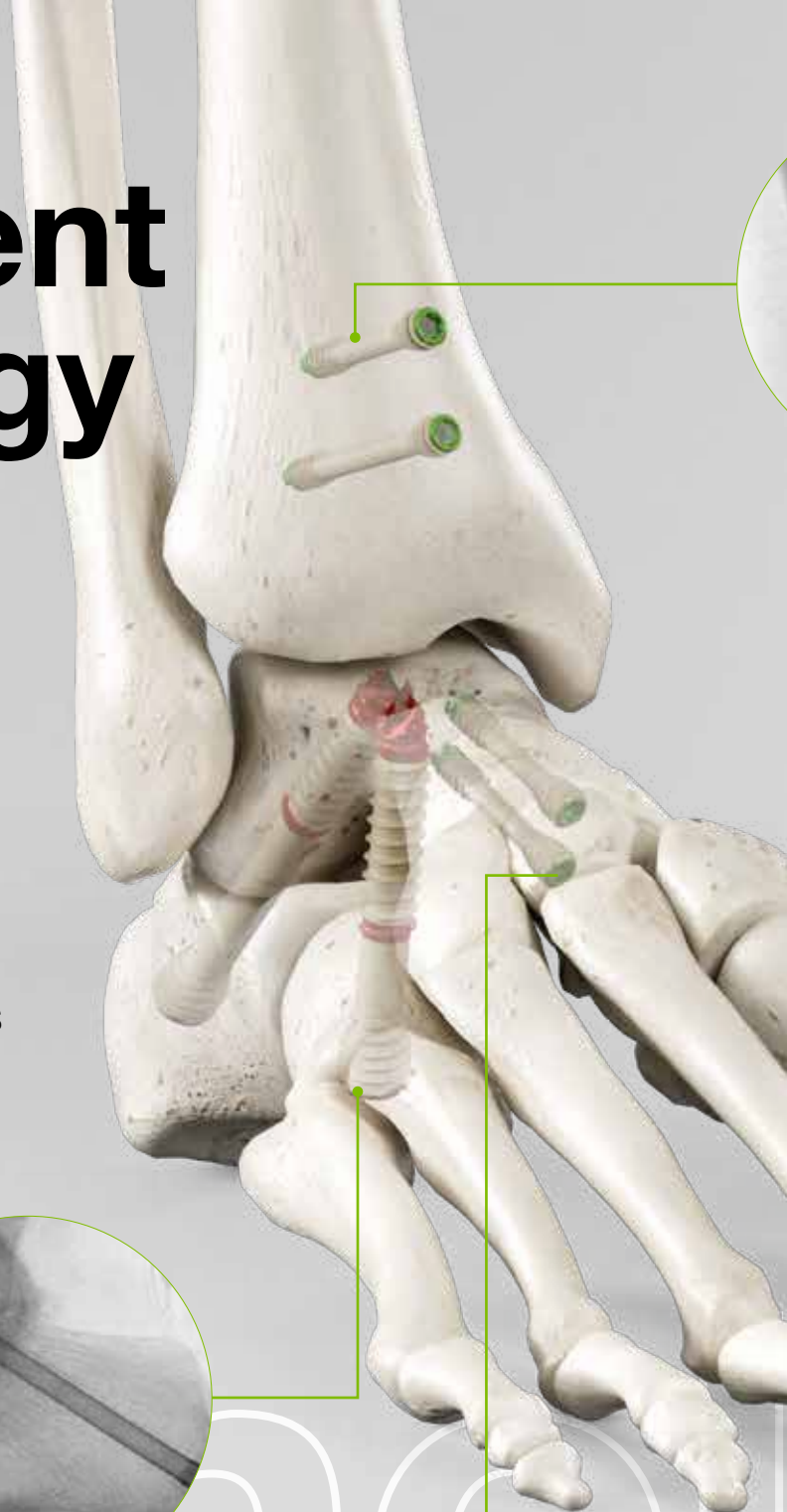

Subtalar fusion

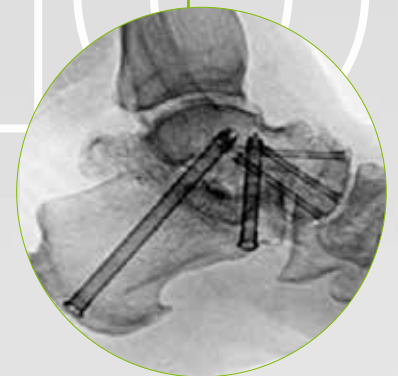

Medial double arthrodesis (TN and Subtalar fusion)

³Patent US11253304B2.


Ortholucen




Posterior malleolar fracture



Ankle fusion

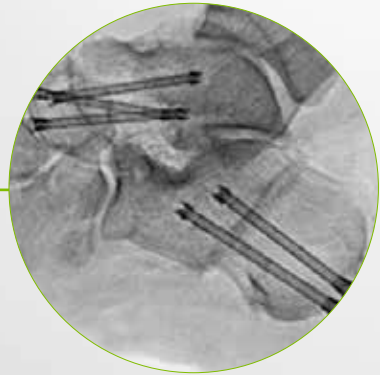

Cavovarus correction
(Dwyer calcaneal osteotomy
and TN fusion)

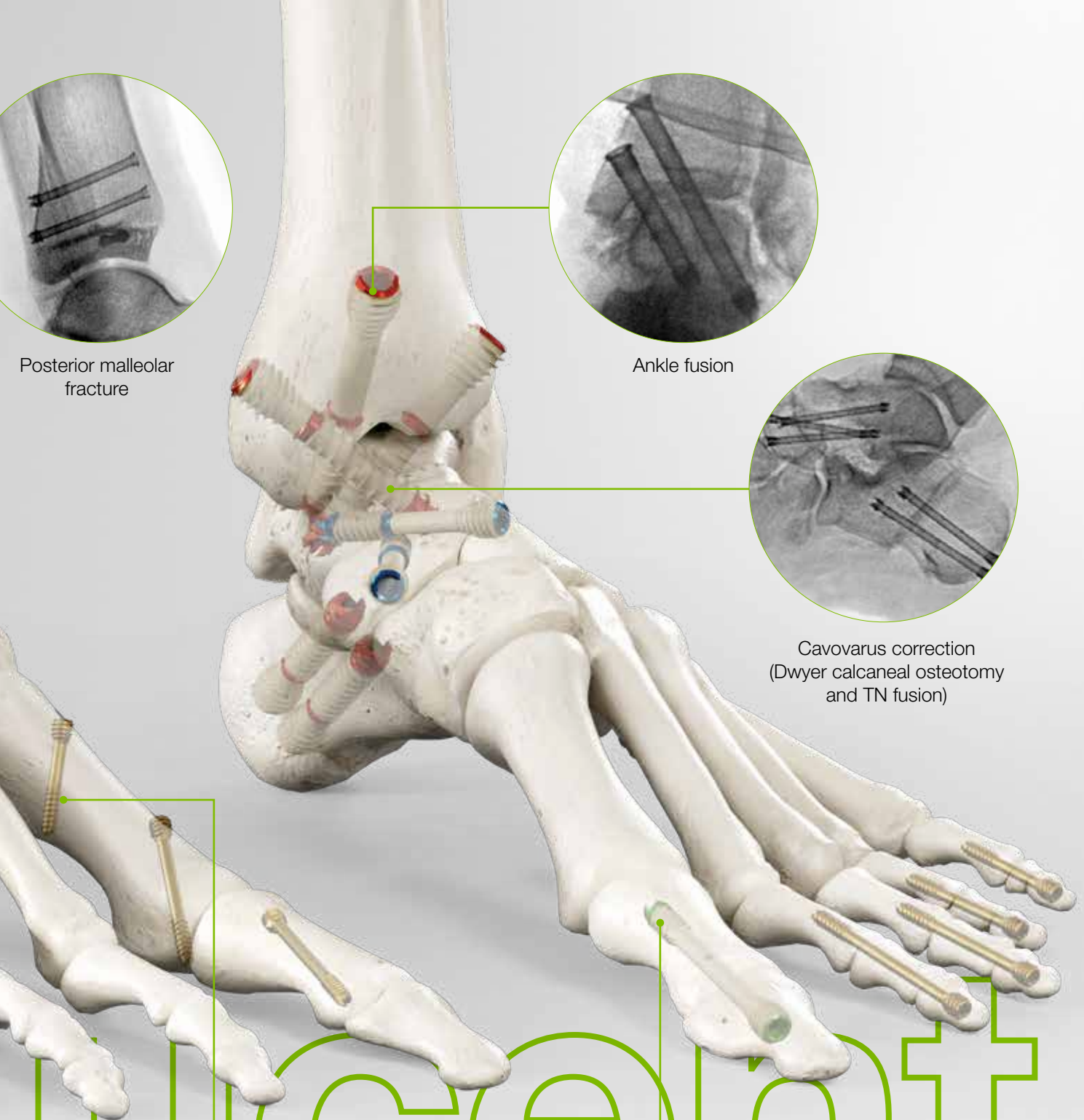

Scarf bunionectomy
with Akin osteotomy




Lesser toe deformities

Wucent

Versatile Cannulated Implant System

Full range of headed and headless implant diameters available in 2.5, 4.3, 5.6 and 7.4mm.

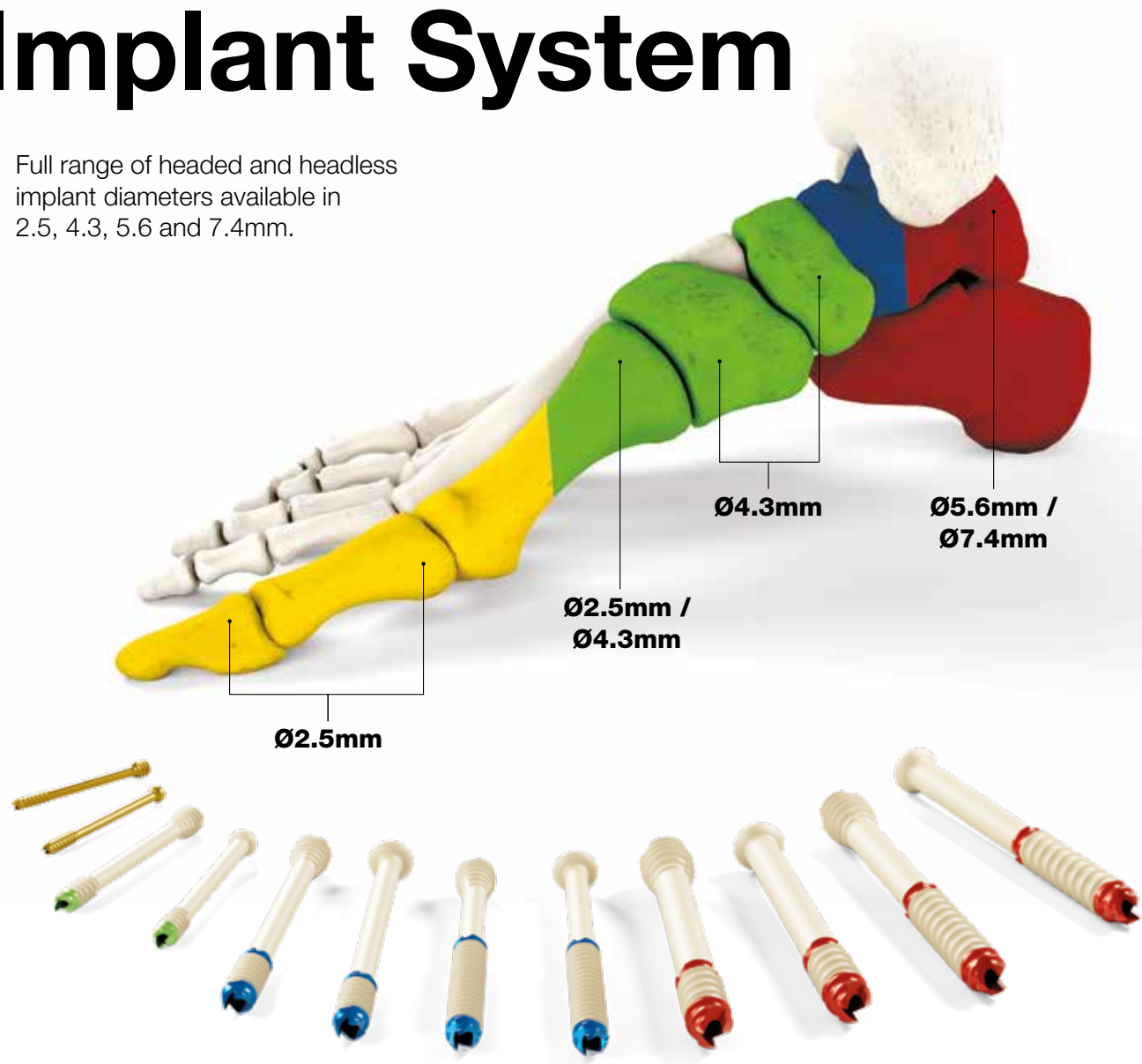

170
Screw
Variations

85
Headless
Screws



85
Headed
Screws

4x
Various
Diameters


5x
Thread
Lengths

Sterile, Disposable Instrument Kits


Single-use

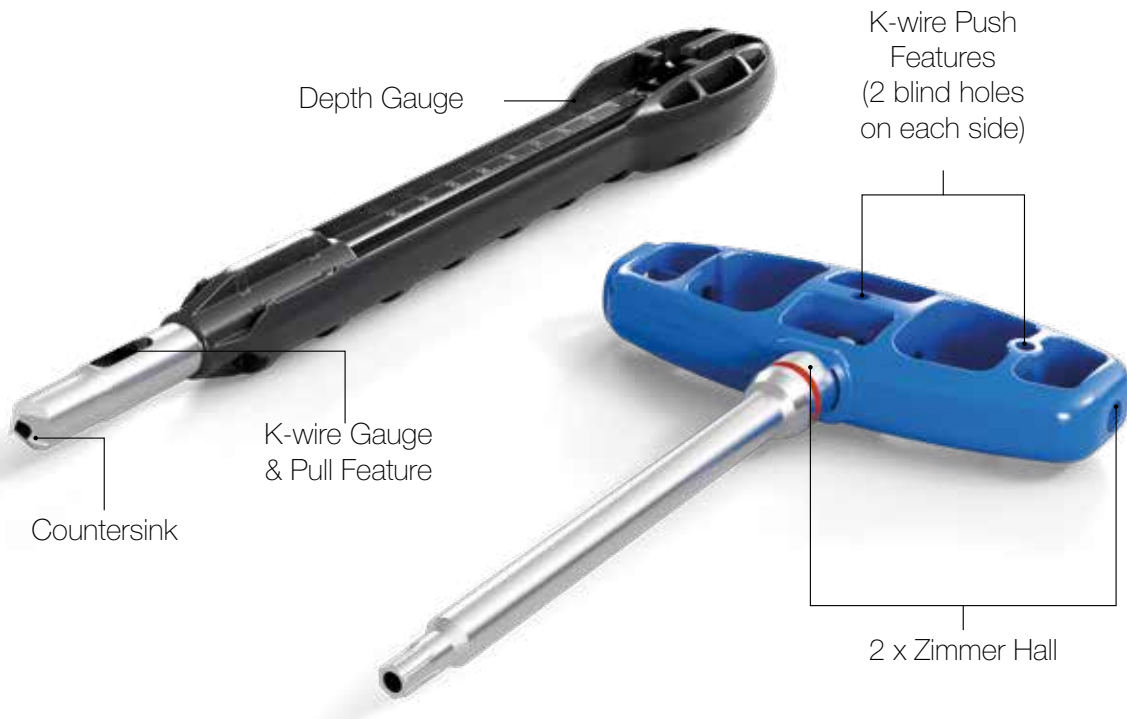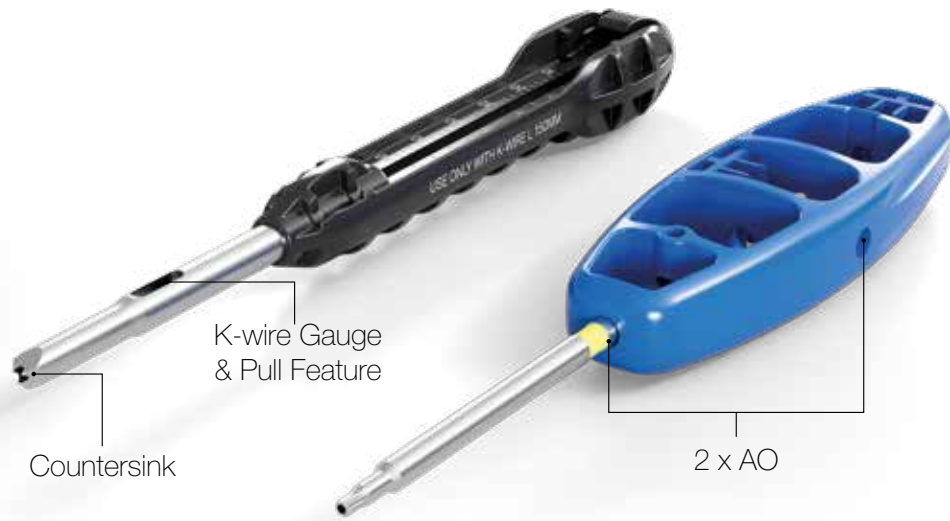

Small Instrument Kit

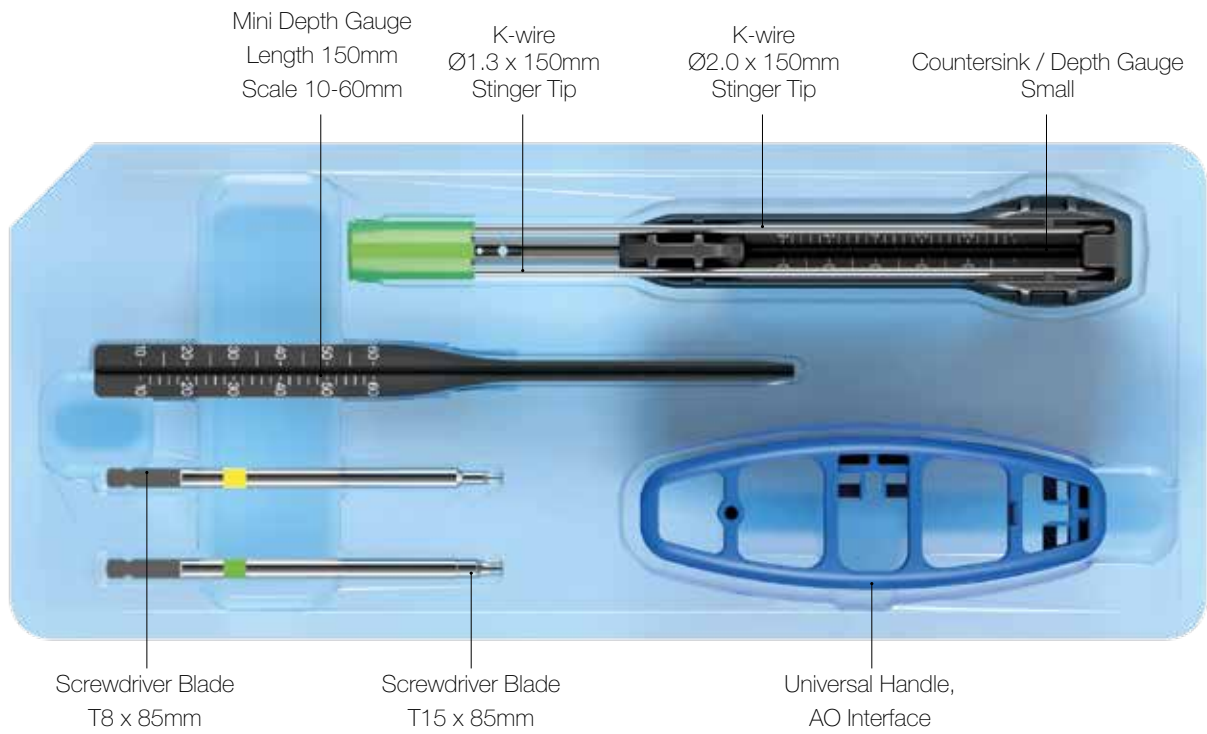

Large Instrument Kit

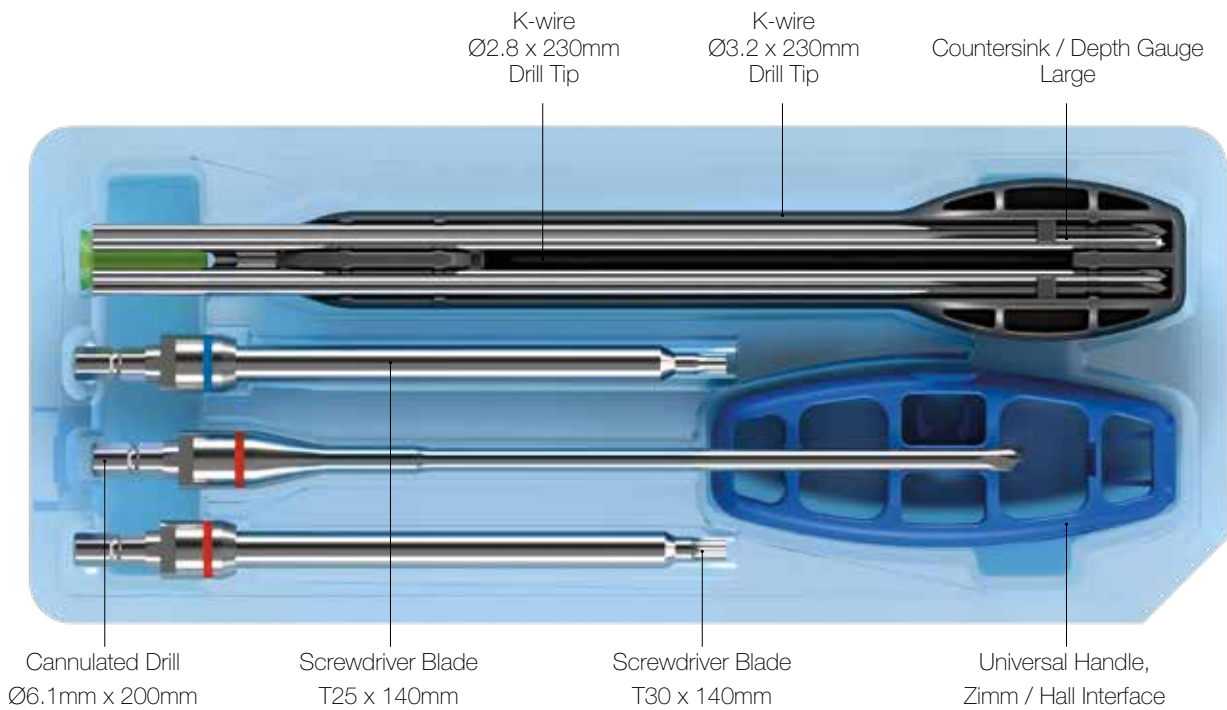

Screw Caddy
with Instrument Kits




Shipping Container



Creed Screws

Instruments Kits

Washers*

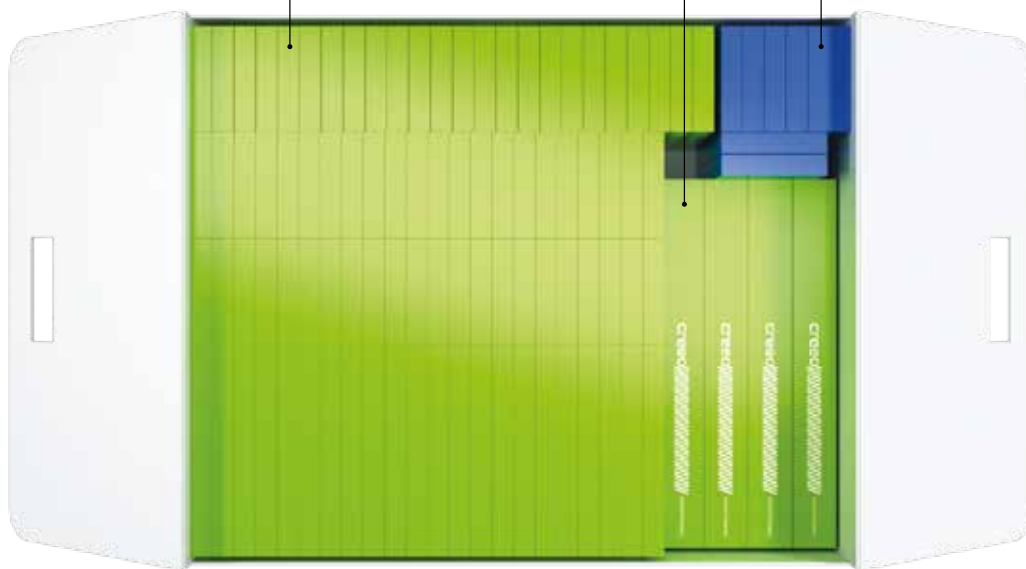

*For headed screw configurations only.



Scan for more
product information

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a surgeon. Rx only.

GLW, Inc, GLW Medical Innovation, and Creed are trademarks of GLW, Inc. Innov8ortho is the exclusive distributor of Creed Ortholucent Implants.

www.glwmed.com

© 2025 GLW, Inc. All rights reserved.
Patents: www.glwmed.com/patents/



Legal Manufacturer:

GLW, Inc.
930 Sylvan Ave, Suite 125
Englewood Cliffs, NJ 07632



Distributed by:

Innov8ortho, LLC
930 Sylvan Ave, Suite 125
Englewood Cliffs, NJ 07632
custsvc@innov8ortho.com

Ref: CRE-BRO-01-25-EN